NEW K-Tape The Original 4 Rolls Brochure FOR PACKAGING PRINT.pdf 1 2020-03-06 1:53 PM

NEW K-Tape The Original 4 Rolls Brochure FOR PACKAGING PRINT.pdf 2 2020-03-06 1:53 PM

NEW K-Tape The Original 4 Rolls Brochure FOR PACKAGING PRINT.pdf 3 2020-03-06 1:53 PM

The Original K-Tape created by Kumbrink is recognized worldwide for its optimum quality and performance. Often copied but never matched.

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
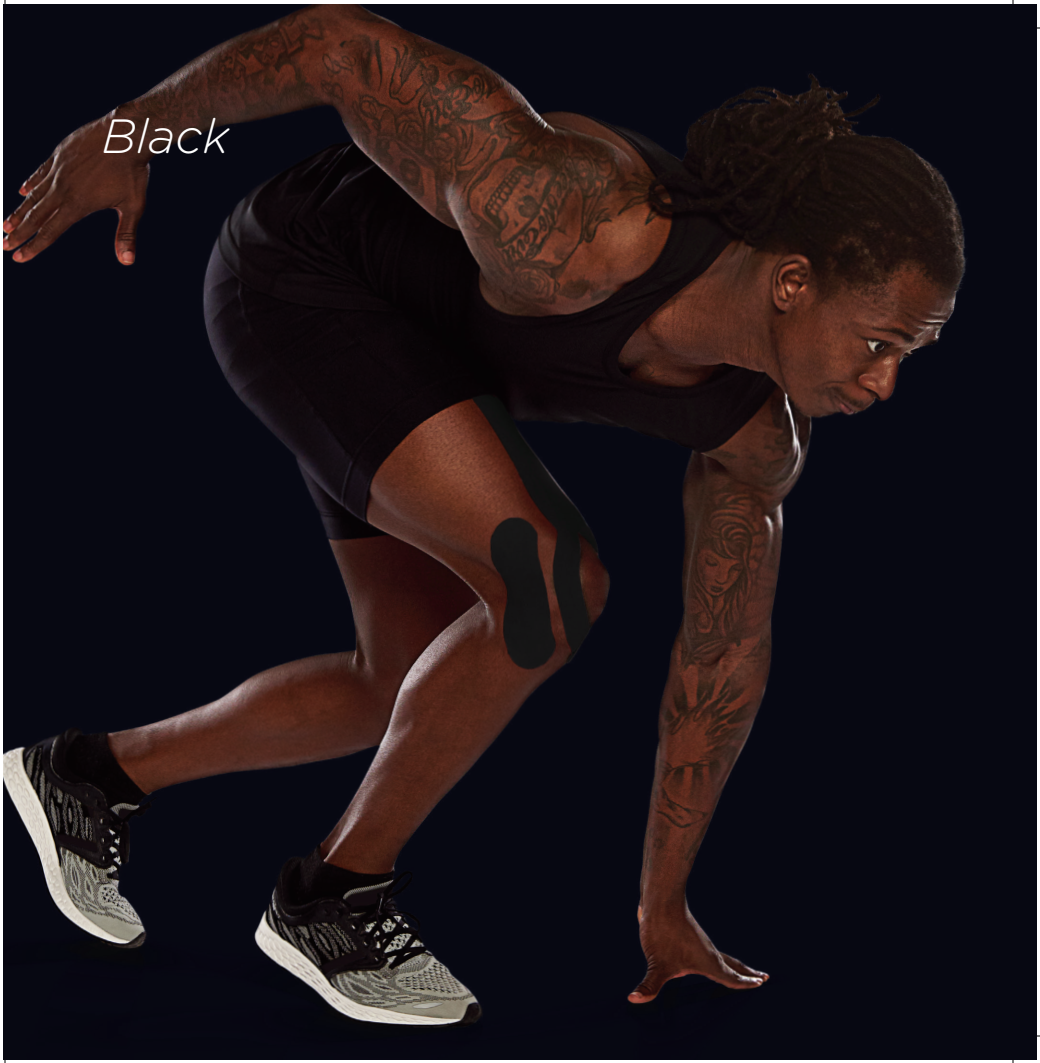
The Original

- *Elastic cotton tape for professional use*
- *Water resistant*
- *Sweat-permeable & breathable*
- *Medication free*
- *No active ingredients*
- *Physiobond acrylic glue*
- *96% cotton, 4% spandex*

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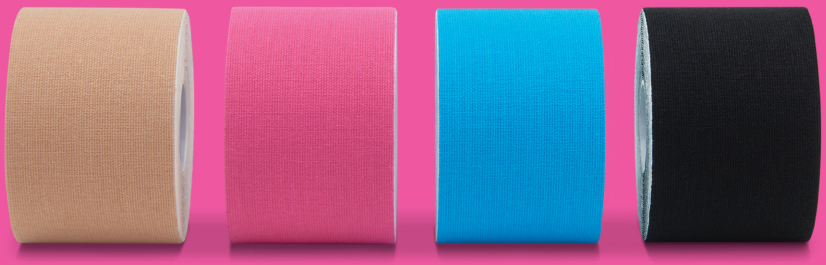K-Tape is an innovative product, recognized worldwide for its optimum quality and performance. The Kumbrink developed original K-Tape combines the best natural, raw materials with the most advanced manufacturing process available.

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Notes

- K-Tape is applied to the skin via a skin-friendly acrylic coating and can be worn for several days.
- K-Tape should only be applied by trained K-Taping therapists.
- In rare cases skin irritation or discoloration can occur.
- Improper use can cause skin lesions.
- If the skin does react to the tape, remove the tape immediately and consult a physician if necessary.
- Apply 30-60 minutes before sporting activities.
- Do not apply to skin affected by any skin condition or irritation.
- Do not apply to open wounds or bleeding areas of the body.
- Do not dry the tape with a hairdryer following showering or bathing.
- Swab wet tape dry. Do not rub as this may loosen tape ends.
- Non-sterile.
- For single use only.

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© K-Tape, K-Tape My Skin and Physiobond are registered trademarks in Germany, Europe, USA, Canada and many more countries. All rights reserved.
biviax GmbH & Co. KG