
KUMBRINK

CROSSTAPE®

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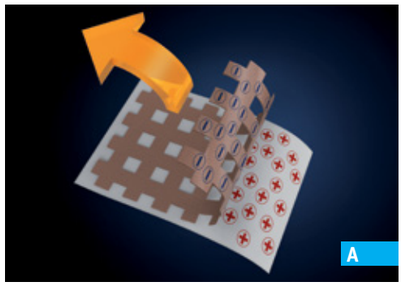KUMBRINK-CROSSTAPE® AND THE CROSSPUNCTURE® THERAPY EN

KUMBRINK-CROSSTAPE® can be used for local points of pain, trigger points, tense muscles, painful joints, headaches, painful scars and many more indications. KUMBRINK-CROSSTAPE® is also the fundamental component of the Crosspuncture® Therapy. KUMBRINK-CROSSTAPE® is used to stimulate selected acupuncture points, enabling an acupuncture-like treatment without the use of needles.

The water-resistant CROSSTAPEs are free of medication and active ingredients. The tapes are applied directly over points of pain, muscular trigger points and acupuncture points. Depending on the stresses applied (such as showering, swimming, sport or work), they can adhere to the skin for a period up to several days.

KUMBRINK-CROSSTAPE® are included in the sizes M, L and XL. The size M is used most of the time. Size L, XL or several size M tapes can be used on larger areas of pain. When doing so, it is important to remember that KUMBRINK-CROSSTAPE® are not elastic and are unable to follow the skin when it is stretched. It is therefore preferable to use several size M KUMBRINK-CROSSTAPE® on areas where the skin is subject to strong stretching. The larger tapes are suitable for parts of the body less subject to stretch, such as the outer surfaces of the thigh and shoulder or the shin, etc.

How does the human body function?

Very small electrical impulses control a large number of functions and information processing in the human body. Whether muscular, fascial or neural functions or acupuncture points, the body frequently uses measurable electrical circuits and resistance. The skin is an important organ here. The largest human organ, equipped with a multitude of receptors and complexly connected to the inner organs and the brain, is able to farther transmit information, pain and perceptive input. Whether for diagnostics (ECG) or treatment (electrostimulation devices), the medical field has already been taking advantage of how the skin functions and its connections for many years now. The same applies to acupuncture, where the exact location of the acupuncture points to be treated are found via their electrical charge. Injuries, diseases, scars and tense muscles — all of these affect the body's electrical conduction system and result in dysfunction and pain being transmitted to the brain.

Properties of KUMBRINK-CROSSTAPE®

KUMBRINK-CROSSTAPE® is made of dual mixed fibres and is attached to a specially coated backing paper. When pulled off the backing paper, the KUMBRINK-CROSSTAPE® are charged up electrostatically so that they are suitable for Crosspuncture® Therapy. This means that the KUMBRINK-CROSSTAPE® has a surplus of electrical charge once pulled off the backing paper that cannot discharge by itself **(A)**. Once the KUMBRINK-CROSSTAPE® has been charged in this manner, it is applied to the skin to stimulate the areas of pain and acupuncture points where it slowly releases its electrical charge.

How to use KUMBRINK-CROSSTAPE®

It is very easy to use KUMBRINK-CROSSTAPE®. Carefully pull the KUMBRINK-CROSSTAPE® off the backing paper **(B)** and let a corner of the tape stick to your fingertip **(C)**. In doing so, touch the KUMBRINK-CROSSTAPE® as little as possible as some of the charge may be lost when the tape is touched repeatedly. When you position the KUMBRINK-CROSSTAPE® over an acupuncture or trigger point approx. one centimetre parallel to the skin, you can often see that the charged tape is attracted to/pulled towards this oppositely charged area **(D)**.

KUMBRINK-CROSSTAPE® can be applied to any point of pain or be used instead of needles over the entire acupuncture spectrum. The exact position of the acupuncture points is important in needle-based acupuncture treatment. As KUMBRINK-CROSSTAPE® cover a larger area, the position does not need to be as accurate, meaning that laypeople are able to apply the tapes to acupuncture points. In the following section you can see examples of use for common indications. Further information can be found on the www.crosstape.com website.

Contraindications

There are no known contraindications or interactions. Nevertheless, we ask that you please pay attention to the following points.Do not apply tape to wounds or skin injuries. Do not use on skin diseases and irritation such as neurodermatitis or psoriasis. Test tolerance on sensitive skin first with a single KUMBRINK-CROSSTAPE®. If the skin does react to the tape, remove the tape immediately and consult a physician if necessary.

Further application techniques and information on KUMBRINK-CROSSTAPE® and Crosspuncture® Therapy can be found at www.crosstape.com

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EN Further application techniques and information on KUMBRINK-CROSSTAPE® and Crosspuncture® Therapy can be found at:

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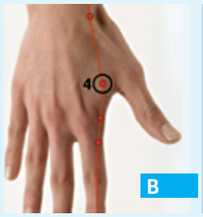
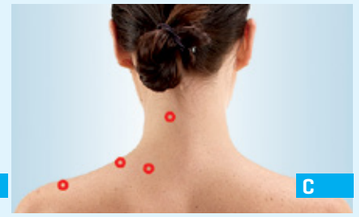
HEADACHE

FOR HEADACHE IN THE TEMPLES
- Apply one size M KUMBRINK-CROSSTAPE® directly onto the left temple and the right temple accordingly **(A)**.
- Apply one KUMBRINK-CROSSTAPE® on the LI-4 acupuncture points on the left and right hands accordingly, between the thumbs and the index fingers **(B)**.

HEADACHES ARISING FROM TENSION IN THE NECK
- Depending on the size of the painful point/area, apply size M or L KUMBRINK-CROSSTAPE® directly over the palpable points of pain in the shoulder/neck region **(C)**.
- Apply one KUMBRINK-CROSSTAPE® on the LI-4 acupuncture points on the left and right hands accordingly, between the thumbs and the index fingers **(B)**.

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
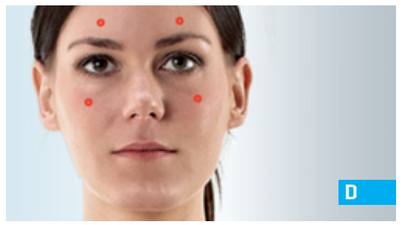
SINUSITIS
- Apply one size M KUMBRINK-CROSSTAPE® over each of the sinuses (centred somewhat over the eyebrows) and one over each of the cheek bones (centrally) **(D)**.
- Apply one KUMBRINK-CROSSTAPE® on the LI-4 acupuncture points on the left and right hands accordingly, between the thumbs and the index fingers **(E)**.

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
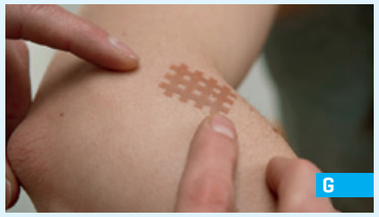
EPICONDYLITIS („TENNIS ELBOW")
- Depending on the size of the painful point/area, apply size M or L KUMBRINK-CROSSTAPE® directly over the palpable point of pain at the elbow **(F)**.
- Apply one size M KUMBRINK-CROSSTAPE® over the LI-4 acupuncture point of the affected arm. LI-4 can be found between the thumb and the index finger **(G)**.

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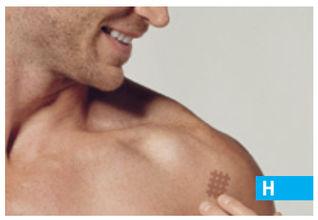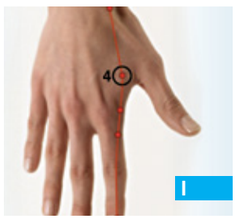SHOULDER PAIN
- Depending on the size of the painful point/area, apply size M, L or XL KUMBRINK-CROSSTAPE® to the point or area of pain at the shoulder **(H)**.
- Apply one size M KUMBRINK-CROSSTAPE® over the LI-4 acupuncture point of the affected arm. LI-4 can be found between the thumb and the index finger **(I)**.
- Additionally apply a size M or L KUMBRINK-CROSSTAPE® so that the Lu-7 and Lu-8 **(J)** acupuncture points are covered over. The exact location can be measured 2 fingers from the edge of the wrist. Tip: The pulse can be easily felt over Lu-7.

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STOMACH PROBLEMS
- Apply one size M KUMBRINK-CROSSTAPE® to each of the Sp-3, Sp-4 and Sp-6 acupuncture points **(K)**. Sp-3 is located directly behind the ball of the big toe, Sp-4 approx. in the middle between the base of the big toe and the start of the heel. Sp-6 is found approx. 3 finger widths above the upper edge of the ankle.
- Apply one size L KUMBRINK-CROSSTAPE® over the CV-15 **(L)** acupuncture point. This is found just underneath the tip of the breast bone and can be directly felt most of the time when the stomach aches.

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