INSTRUCTIONS FOR USE

EMBRACE PRE-FILLED FLUSH SYRINGE

0.9% sodium chloride injection


- PRE-FILLED FLUSH SYRINGE
- VORGEFÜLLTE SPÜLSPRITZE
- SERINGUE DE RINÇAGE PRÉ-REMPLIE
- JERINGA DE IRRIGACIÓN PRE-LLENADA
- SIRINGA DI LAVAGGIO PRERIEMPITA
- Seringa pré-cheia flush
- VOORGEVULDE SPOELSPUIT
- ELŐRETÖLTÖTT ÖBLITŐ FECSKENDŐ
- FÖRFYLLD SPOLSPRUTA
- FORFYLDT SKYLLESPRØJTE
- ESITÄYTETTY HUUHTELURUISKU
- WSTĘPNIE NAPEŁNIONA STRZYKAWKA DO PRZEPŁUKIWANIA
- PŘEDNAPLNĚNÁ PROPLACHOVACÍ STŘÍKAČKA
- PREDNAPLNENÁ PREPLACHOVACIA STRIEKAČKA
- ПРЕДВАРИТЕЛНО НАПЪЛНЕНА СПРИНЦОВКА ЗА ПРОМИВАНЕ
- PREDNAPOLNJENA BRIZGA ZA IZPIRANJE
- ÖNCEDEN DOLDURULMUŞ YIKAMA ENJEKTÖRÜ
- ПОПЕРЕДНЬО НАПОВНЕНИЙ ШПРИЦ ДЛЯ ПРОМИВАННЯ
- ΠΡΟΓΕΜΙΣΜΕΝΗ ΣΥΡΙΓΓΑ ΠΛΥΣΗΣ
- ÖNCEDEN DOLDURULMUŞ YIKAMA ENJEKTÖRÜ

1

PRE-FILLED FLUSH SYRINGE

PRODUCT STRUCTURE:

1. Protective Tip Cap - Luer Lock
2. 0.9% sodium chloride injection
3. Syringe Barrel
4. Piston
5. Plunger Rod

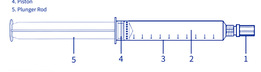INTENDED USE:
The device is used for closing and flushing the end of catheters at intervals of different drug treatment.

CONTRAINDICATIONS: None

PRECAUTIONS:

1. The use of the product must meet the requirements of relevant operational specifications and relevant regulations of medical department, and be restricted to trained doctors or nurses.
2. Check for package integrity. It is strictly prohibited to use if product package is damaged, or the syringe breaks, leaks, disintegrates.
3. It is strictly prohibited to use if protective tip cap is not installed correctly, loose or falls off.
4. Visual examination for clarity of solution, if solution discolouration, turbidity, precipitation or any form of suspended particulate matter material appears-DO NOT USE.
5. DO NOT RE-STERILIZE
6. Please use the product within the validity period. The product batch number, production date and expiration date are shown on the product packaging label.
7. It is strictly prohibited to re-use. Single use only, and discard any unused remaining solution.
8. Check the compatibility of 0.9% sodium chloride solution with drugs. Direct contact with drugs incompatible with 0.9% sodium chloride injection is strictly prohibited.
9. It is strictly prohibited to be used for subcutaneous and intramuscular injection and dispensing. Only for flushing of IV-set catheters. It is not for medication dilutions, nor for dry product reconstitution.
10. DO NOT PLACE THE PRODUCT ON A STERILE FIELD.
11. Do not use if you are allergic to the sodium chloride solution.
12. Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
13. If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.
14. Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

⚠	Caution	⊘	Do not use if package is damaged
EC REP	Authorized representative in the European Community	LOT	Lot number
CE 0123	CE Mark: conforms to essential requirements of the Medical Device Directive 93/42/EEC	⊘	Do not re-sterilize
	Date of manufacture		Consult instructions for use
	Manufacturer		Keep away from sunlight
	Use-by date	REF	Catalogue number
	Keep dry	STERILE R	Sterile irradiation using ionizing beams
	Do not re-use		Store at controlled temperature(20-30°C).

Jiangsu Embrace Science & Technology Development Co., Ltd.
29-2 Yinghin Road, Huai'an Economic and Technological Development Zone, Jiangsu, Zip Code: 223001
www.emtdi.cn

EC REP **Riomavis S.L.**
Calle de Almansa 55, 1D, Madrid 28039 Spain
Tel.: +34 658 396 230

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CH REP **MedNet SWISS GmbH**